beurer

cellulite releaZer®

DE	cellulite releaZer compact Gebrauchsanweisung	2
EN	cellulite releaZer compact Instructions for use	21
FR	cellulite releaZer compact Mode d'emploi	40
ES	cellulite releaZer compact Instrucciones de uso	59
IT	cellulite releaZer compact Istruzioni per l'uso	78
TR	cellulite releaZer compact Kullanım kılavuzu	97
RU	cellulite releaZer compact Инструкция по применению	116
PL	cellulite releaZer compact Instrukcja obsługi	136
NL	cellulite releaZer compact Gebruiksaanwijzing	155
DA	cellulite releaZer compact Betjeningsvejledning	174
SV	cellulite releaZer compact Bruksanvisning	192
NO	cellulite releaZer compact Bruksanvisning	211
FI	cellulite releaZer compact Käyttöohje	229

DEUTSCH

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ENGLISH

Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

WARNING

- The device is only intended for domestic/private use, not for commercial use.
- This device may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the consequent risks of use.
- Children must not play with the device.

- If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of.
- Under no circumstances should you open or repair the device yourself, as faultless functionality can no longer be guaranteed thereafter. Failure to adhere to this instruction will void the warranty.
- Cleaning and user maintenance must not be performed by children unless supervised.

Table of contents

1. Included in delivery	24
2. Signs and symbols	24
3. Notes	26
4. Intended use	28
5. Device description	30
6. Initial use/charging process	31
7. Use.....	32
7.1 Preparing for use	32
7.2 Switching on the cellulite releaZer compact / selecting vibration setting	33
7.3 Tips for use.....	34
7.4 Treatment on the outside of the thigh.....	34
7.5 Treatment on the rear of the thigh	36
8. Cleaning and maintenance.....	37
9. Disposal	38
10. Technical data	39
11. Warranty / Service	39

1. Included in delivery

Check that the device packaging has not been tampered with and make sure that all contents are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Services address.

1x cellulite releaZer compact / 1x micro-USB charging cable /
1x mains adapter / 1x these instructions for use

2. Signs and symbols

The following symbols appear in these instructions for use:

	WARNING Warning notice indicating a risk of injury or damage to health
	IMPORTANT Safety note indicating possible damage to the device/accessory
	Note Note on important information
	Observe the instructions for use

	Manufacturer
	Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive – WEEE
	Separate the packaging elements and dispose of them in accordance with local regulations
	Marking to identify the packaging material A = Material code, B = Material number: 1-7 = Plastics, 20-22 = Paper and cardboard
	Only use indoors
	The device is double protected and therefore corresponds to protection class 2
IP20	Protected against solid foreign objects 12.5 mm in diameter and larger
	Energy efficiency class 6
	Direct current
	Short-circuit-proof safety isolating transformer

	Switching mains part
	Polarity.
	CE labelling This product satisfies the requirements of the applicable European and national directives.
	The products demonstrably meet the requirements of the Technical Regulations of the EAEU.
	United Kingdom Conformity Assessed Mark
	Separate the product and packaging elements and dispose of them in accordance with local regulations.

3. Notes

WARNING

Safety notes regarding the use of the cellulite releaZer compact

- Keep the packaging material away from children. Choking hazard!
- The cellulite releaZer compact is not a medical device. It is a massage treatment device.

- Before use, ensure that the device and accessories have no visible damage. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.
- Repairs must only be carried out by Customer Services or authorised retailers. Under no circumstances should you open or repair the cellulite releaZer compact yourself, as faultless functionality can no longer be guaranteed thereafter. Failure to comply with this instruction will void the warranty. Never attempt to repair the device yourself!
- Only use the device with the supplied accessories.

Notes on handling batteries

- If your skin or eyes come into contact with fluid from a rechargeable battery cell, flush out the affected areas with water and seek medical assistance.
- **⚠ Choking hazard!** Small children may swallow and choke on batteries. Store batteries out of the reach of small children.
- Protect batteries from excessive heat.
- **⚠ Risk of explosion!** Never throw batteries into a fire.
- Do not disassemble, split or crush the batteries.
- Only use chargers specified in the instructions for use.
- Batteries must be charged correctly prior to use. The instructions from the manufacturer and the specifications in these instructions for use regarding correct charging must be observed at all times.
- Fully charge the battery prior to initial use (see Chapter 7. Initial use/charging process).
- In order to achieve as long a battery service life as possible, fully charge the battery at least twice per year.

4. Intended use

WARNING

Only use the cellulite releaZer compact for massage treatment of cellulite. The cellulite releaZer compact is intended for private use only. The cellulite releaZer compact is not a medical device; rather it is a massage treatment device. The cellulite releaZer compact may only be used for the purpose for which it was designed and in the manner specified in the instructions for use. Any form of improper use can be dangerous. The manufacturer is not liable for damage resulting from improper or incorrect use.

- Do not use the cellulite releaZer compact if you suffer from a medical abnormality or injury in the area to be massaged (e.g. open wounds).
- Do not use the cellulite releaZer compact during pregnancy.
- Do not use the cellulite releaZer compact whilst sleeping.
- Do not use the cellulite releaZer compact in a vehicle.
- Do not use the cellulite releaZer compact on animals.
- Do not use the cellulite releaZer compact whilst undertaking any activity where an unexpected reaction could be dangerous.
- Do not use the cellulite releaZer compact on your head.
- Do not use the cellulite releaZer compact on your face.
- Do not use the cellulite releaZer compact on your genitals.
- Do not use the cellulite releaZer compact after the consumption of substances that lead to impaired cognitive function (e.g. pain medication, alcohol).
- Do not use the cellulite releaZer compact for longer than 15 minutes (risk of overheating) and leave the device to cool down for at least 15 minutes before using it again. The cellulite releaZer compact has an automatic switch-off function meaning that it switches off automatically after 15 minutes of use.


- Do not use the cellulite releaZer compact either before or after sunbathing. Doing so increases the risk of sunburn.

If you are unsure whether the cellulite releaZer compact is suitable for you, please consult your doctor.

Please consult your doctor before using the cellulite releaZer compact

- if you suffer from a severe illness or have recently had surgery,
- if you have a pacemaker, implant or other aids,
- in the case of feverish colds, varicose veins, thrombosis, phlebitis, jaundice, diabetes, nerve disorders (e.g. sciatica) or acute inflammations,
- if you are from suffering pain, the cause of which is unknown.

5. Device description

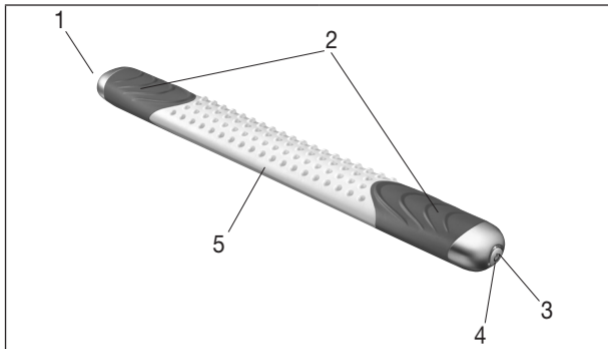1. Micro USB connection

2. Handles

3. ON/OFF button

4. Function light

5. Four treatment surfaces
(smooth, nubs, blunt edge, and sharp edge)

Function light status 	Charging function (charging cable plugged in)
Function lights flash in a circle	Charging
Function lights illuminate continuously	Charging complete
The left and right function lights flash alternately	Low battery

6. Initial use/charging process

Note

Before using the cellulite releaZer compact for the first time, let it charge for a minimum of 3.5 hours. Proceed as follows:

- Connect the charging cable to the mains adapter and the cellulite releaZer compact (*Figure 1*). Then insert the mains adapter into a suitable socket. Alternatively, you can also charge the cellulite releaZer compact with the charging cable via the USB port on your computer or laptop (*Figure 2*).

When the device is charging, the function lights flash in a circle. When the cellulite releaZer compact is fully charged, the function lights illuminate continuously. If the left and right function lights flash alternately during operation, you have to charge the cellulite releaZer compact.

Figure 1


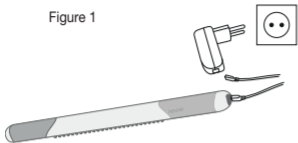
Figure 2

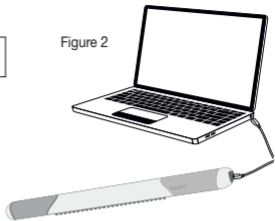7. Use

7.1 Preparing for use

Before you start using the device, you should gently rub massage oil onto the parts of the body to be treated.

Do not use any oils or creams with a strong cooling or heating effect, as they may cause skin irritation.

1. Put a couple of drops of massage oil onto the palm of your hand.
2. Now rub the massage oil evenly onto the parts of the body to be massaged.
3. Before treatment, wipe your hands with a towel to prevent them from slipping when using the device and to prevent damage to the silicone surface on the handles.

7.2 Switching on the cellulite releaZer compact / selecting vibration setting

Press the ON/OFF button to switch on the cellulite releaZer compact. The function light now illuminates blue and the first vibration setting is activated. The cellulite releaZer compact has a total of four different vibration settings.

Action	Vibration setting
1st press of the button	The cellulite releaZer compact switches on and the function light illuminates blue. Low vibration setting.
2nd press of the button	The second function light illuminates blue. Medium vibration setting.
3rd press of the button	The third function light illuminates blue. High vibration setting.
4th press of the button	The fourth function light flashes blue. Intermittent vibration setting.
5th press of the button	The cellulite releaZer compact switches off and the function light goes out.

7.3 Tips for use

- During use, hold the cellulite releaZer compact with both hands using the handles.
- For optimum results, use the cellulite releaZer compact every two days.
- Severe reddening or bruises may occur at the treated areas after use. This is completely normal and indicates that you have used the cellulite releaZer compact correctly.

Note


You can use the cellulite releaZer compact on all areas where you have cellulite. The examples in the following chapter showing the outside of the thigh (iliotibial band) and the rear of the thigh (biceps femoris muscle) demonstrate how the cellulite releaZer compact should be used.

7.4 Treatment on the outside of the thigh

The treatment time for the outside of the thigh is approx. 7 minutes and consists of a total of three treatment phases. To start treatment, proceed as follows:

1. Follow the instructions in the "Preparing for use" chapter.
2. Place the leg to be treated on a chair, stool or the edge of the bath tub, for example. Carry out the entire treatment on the thigh whilst standing.
3. Select your desired vibration setting (1, 2, 3 or 4) using the ON/OFF button. The cellulite releaZer compact starts to vibrate.


4. Treatment phase I: Treat the outside of the thigh with the nub side using small, fast, circular movements and moderate pressure (for approx. 3 minutes).
5. Treatment phase II: Treat the outside of the thigh with the nub side using slow, deep, circular movements and increased pressure (for approx. 3 minutes).
6. Treatment phase III: The third treatment phase is drainage. Glide the device towards the centre of the body/heart using the coloured edge and moderate pressure. When doing so, start the movement at the upper part of the outside of the thigh and slowly move the starting point down towards the knee. One stroke should be approximately 10 cm long. Carry out this process 5 times moving down the outside of the thigh until you reach the knee area. At the end, move the device over the tissue of the entire outside of the thigh once.

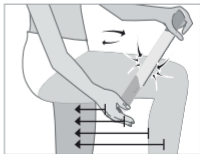
7.5 Treatment on the rear of the thigh

The treatment time for the rear of the thigh and buttocks is approx. 7 minutes and consists of a total of three treatment phases. To start treatment, proceed as follows:

1. Follow the instructions in the "Preparing for use" chapter.
2. Place the leg to be treated on a chair or stool. Carry out the entire treatment on the thigh whilst standing.
3. Select your desired vibration setting (1, 2, 3 or 4) using the ON/OFF button. The cellulite releaZer compact starts to vibrate.
4. Treatment phase I: Treat the rear of the thigh with the nub side using small, fast, circular movements and moderate pressure (for approx. 3 minutes).
5. Treatment phase II: Treat the rear of the thigh with the nub side using slow, deep, circular movements and increased pressure (for approx. 3 minutes).
6. Treatment phase III: The third treatment phase is drainage. Glide the cellulite releaZer compact towards the centre of the body/heart using the coloured edge and moderate pressure. When doing so, start the movement at the upper part of the rear of the thigh (buttocks) and slowly move the starting point down towards the knee.

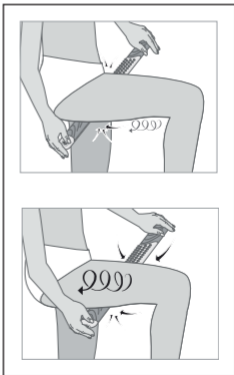The first stroke should be approximately 10 cm long. Carry out this process 5 times moving down the rear of the thigh until you reach the knee area, extending each successive stroke by a couple of centimetres. At the end, move the device over the tissue of the entire rear of the thigh once.

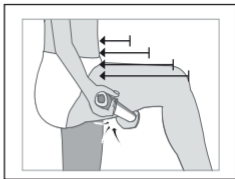8. Cleaning and maintenance

The service life of the cellulite releaZer compact depends on careful use.

IMPORTANT

- Clean the cellulite releaZer compact with a wet or dry cloth after each use. Do not use any abrasive cleaning products.
- Protect the cellulite releaZer compact from knocks, dust, chemicals, drastic changes in temperature, electromagnetic fields and nearby sources of heat (ovens, heaters).

9. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life.

Dispose of the device at a suitable local collection or recycling point in your country. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

Battery disposal

The empty, completely flat batteries must be disposed of through specially designated collection boxes, recycling points or electronics retailers. You are legally required to dispose of the batteries. To remove the battery from the device, pry open the cover on the micro USB connection using a screwdriver. Loosen the four screws on the battery holder. Then remove the battery.

The codes below are printed on batteries containing harmful substances:

Pb = Battery contains lead

Cd = Battery contains cadmium

Hg = Battery contains mercury

10. Technical data

Supply Input: Output:	100-240 V~, 50/60 Hz, 0.2 A 5.0 V === 1.0 A, 5.0 W
Average active efficiency	≥ 73.88 %
No-load power consumption	≤ 0.09 W
Dimensions	40 x 3.2 x 4.7 cm
Weight	Approx. 330g
Battery: Capacity Nominal voltage Type designation	2200 mAh 3.7 V Lithium-ion

11. Warranty / Service

Further information on the warranty and warranty conditions can be found in the warranty leaflet supplied.

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BEURER GmbH • Söflinger Str. 218 • 89077 Ulm (Germany) • www.beurer.com